


Open Loop Display Behavior (valid for FreeAPS X 0.2.0)


- It is important to use open loop from time to time to check settings
- Make sure your Nightscout site is connected and that your profile is correct
 - To modify a profile on NS: Menu -> Profile Editor -> Select Duplicate Icon (circled in yellow in graphic) to copy the current profile and then adjust default date/time if necessary and edit settings to match

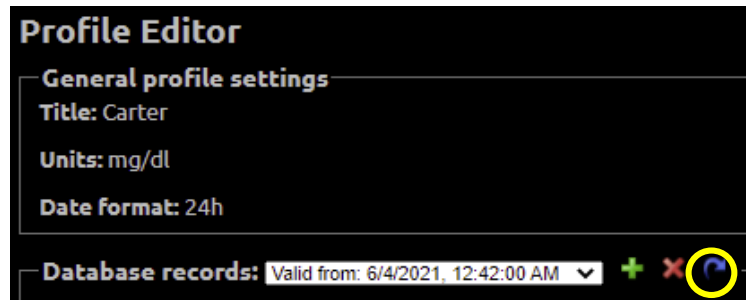

- When you are in open loop, some display features are interesting (HUD is on Phone, NS is on Nightscout – ensure your profile is up to date)
 - If you restore scheduled basal rate, HUD does not show rate but NS does
 - If you enact manual basal rates (explicit value), both HUD and NS plot this
 - When manual TB expires, HUD does not show rate but NS does
 - IOB decays normally, both HUD and NS show this in current pill
 - OpenAPS pill does not include current IOB, so can't see it in history
 - Grabbed some screenshots to show how displays change while in open loop

HUD 05:56

Open, Suggested
Sch basal not plotted

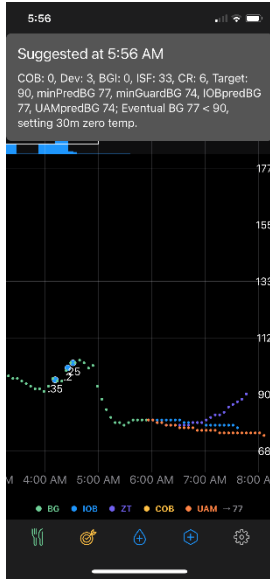


NS 05:57

“Suggested” not
obvious



HUD 06:01


Prev sch basal plots
after TB set

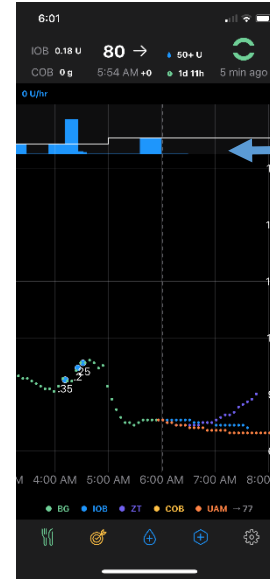

Manual 0 TB
displays on HUB
as tiny blue line

NS 06:32

Immediately shows
scheduled basal
restored when TB
expires.

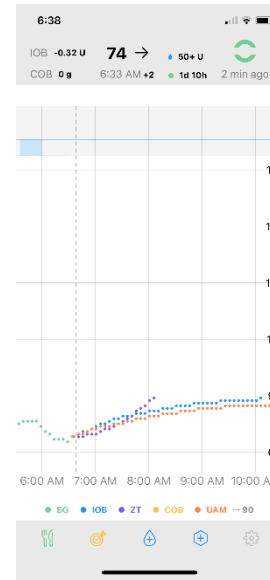

HUD 06:38

Once 0 TB expires, no
longer see tiny blue
line but sch basal is
not plotted